| | | |
|---|---------|---|
| [D] Hautklar-Stift Gebrauchsanweisung | 2 – 5 | [RUS] Прибор для чистки кожи лица Инструкция по применению 26 – 30 |
| [GB] Clear skin pen Instructions for use..... | 6 – 9 | [PL] Sztyft oczyszczający Instrukcja obsługi 31 – 34 |
| [F] Stylo correcteur contre les imperfections Mode d'emploi | 10 – 13 | [NL] Reinigingsstick Gebruiksaanwijzing 35 – 38 |
| [E] Lápiz limpiacutis Instrucciones para el uso..... | 14 – 17 | [P] Caneta limpa-pele Instruções de utilização 39 – 42 |
| [I] Penna depurativa Istruzioni per l'uso | 18 – 21 | [GR] Στικ καθαρισμού προσώπου Οδηγίες χρήσης 43 – 46 |
| [TR] Cilt Temizleme Kalemi Kullanım kılavuzu | 22 – 25 | |

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Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty, baby and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

Included in delivery

- 1 FCE 75 Clear skin pen
- 1 AA 1.5 V battery
- 1 Instructions for use

Intended use

The FCE 75 clear skin pen is to be used exclusively for the treatment of spots and impurities in the human facial area.

Do NOT use the device for severe, nodular or cystic acne.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

This device is intended for domestic use only, not for commercial purposes.

Safety notes**Symbols**

The following symbols appear in these instructions for use.

| | |
|---|---|
|  | Warning instruction indicating a risk of injury or damage to health. |
|  | Note on important information. |
|  | Disposal in accordance with EC Directive 2002/96/EC – WEEE (Waste Electrical and Electronic Equipment). |

Important safety notes **Warning**

Before using the clear skin pen, read the safety notes through carefully. Non-observance of the following information may result in personal injury or material damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the device to third parties.

- Do NOT use the device on moles, warts, ingrown hairs or irritated skin.
- Do not use the clear skin pen if you are particularly sensitive to light.
- Do not use the clear skin pen if you have had sunburn within the past 2 weeks.

- Before use, perform a light compatibility test. To do this, first test the clear skin pen in a less sensitive area such as the inner arm before beginning treatment on the face. If you are unsure, consult your doctor before application.
- Do not apply the treatment to a spot more than 3 times daily. Any more often than this could lead to skin burns or injuries. Leave a minimum time of 4 hours between each treatment.
- If after using the pen you experience a worsening of the skin, undesired side effects or prolonged redness, stop using the clear skin pen and consult your doctor. However, a temporary, mild reddening after application is normal.
- Do not use the clear skin pen in and around the eye and lip area. These areas are particularly sensitive and are therefore not suitable for treatment.
- To prevent damage to the eyes never look directly at the blue light.
- Keep the clear skin pen away from water.
- Do not drop the clear skin pen. The clear skin pen is a precision device. Dropping it from a height could impair the performance of the device.
- The clear skin pen is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised or instructed in the use of the clear skin pen by a person who is responsible for their safety.
- Do not use the clear skin pen if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- Children must not play with the device. Cleaning and user maintenance must not be performed by children unless supervised.

1. Getting to know your device


The FCE 75 clear skin pen is ideal for treating spots and impurities on the face thanks to its blue light therapy.

The clear skin pen is gentle on the skin but effective against impurities. The clear skin pen takes effect within 24 hours after application, enabling spots to subside rapidly. The device is also fitted with a heat function, which promotes circulation.

The practical, pocket size device is perfect for applications at home or on the move.

2. Unit description

1. Treatment surface
2. ON/OFF button
3. Start button
4. Operating indicator (green)
5. Readiness indicator (blue)

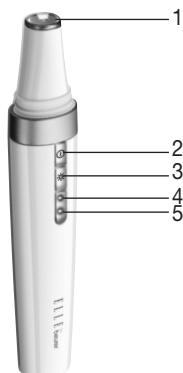3. Operation

- (i)** Do not apply the treatment to a spot more than 3 times daily. You will normally notice an improvement after the first treatment. If necessary, apply the treatment again 4 to 12 hours later.
1. Remove the lid from the clear skin pen.
 - (i)** Clean the heating plate on the clear skin pen carefully before each application using a cloth dampened with soapy water (or alcohol).
 2. Slide the battery compartment lid open. Place the supplied battery in the battery compartment according to the markings on the device. Ensure that the correct battery polarity is observed. Close the battery compartment again.
 3. Press the ON/OFF button.
You will hear a short acoustic signal. The green operating indicator lights up and the blue readiness indicator starts to flash.
 4. After approx. one minute you will hear a number of rapid acoustic signals and the blue readiness indicator will start to flash more quickly. The clear skin pen is now ready for use. To begin the application, press the start button within 15 seconds.
The blue readiness indicator and the blue light on the treatment surface will light up. Place the small heating plate gently on the spot. During the treatment cycle you will hear an acoustic signal every 6 seconds.
- (i)** Do not press the treatment surface down firmly on the spot, as this can lead to severe reddening.

After approx. 2.5 minutes the treatment is complete and you will hear several acoustic signals one after the other. Remove the treatment surface from the affected area.

5. If you want to treat another spot, press the start button again and hold the heating plate to the affected area, otherwise the device will switch off automatically after 15 seconds.

If you do not want to continue the treatment and would like to switch off the device, press the ON/OFF button.

- (i)** A temporary, mild reddening after application is normal.

Warning

- Keep batteries out of the reach of children. Children could put batteries in their mouth and swallow them. This may severely damage their health. In such cases, seek medical assistance immediately!
- Normal batteries must not be charged, heated or thrown into an open fire (risk of explosion!).
- Leaking batteries may damage the device. If you do not intend to use the device for longer periods, remove the battery from the battery compartment.
- Do not use rechargeable batteries!

4. Maintenance and cleaning

- When you are not using the clear skin pen, store it in a cool, dry place away from direct light.
- If you do not intend to use the clear skin pen for longer periods, remove the battery from the battery compartment.
- Clean the heating plate on the clear skin pen carefully using a cloth dampened with soapy water (or alcohol).
- Do not use any cleaning agents or solvents, as they could damage the device.
- Under no circumstances should you hold the device under water, as this can cause liquid to enter and damage the device.

5. What if there are problems?

| Problem | Solution |
|---------------------------|--------------------|
| Device does not switch on | Change the battery |

6. Disposal

Observe the local regulations for material disposal. Dispose of the unit in accordance with EC Directive 2002/96/EC – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

7. Technical data

| | |
|--------------|------------------------|
| Model no. | FCE 75 |
| Power supply | 1 x 1.5 V AA battery |
| Weight | 36 g (without battery) |
| Wavelength | 415 nm, +/-5 nm |

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