|   |       |
|---|-------|
| <b>(D) Ultraschallgerät</b><br>Gebrauchsanleitung                         | 2-11  |
| <b>(GB) Ultrasound device</b><br>Instructions for use                     | 12-21 |
| <b>(F) Appareil à ultrasons</b><br>Mode d'emploi                          | 22-31 |
| <b>(E) Aparato de ultrasonidos</b><br>Instrucciones de uso                | 32-41 |
| <b>(I) Apparecchio a ultrasuoni</b><br>Istruzioni per l'uso               | 42-51 |
| <b>(TR) Ultrason cihazı</b><br>Kullanım kılavuzu                          | 52-60 |
| <b>(RUS) Прибор для ухода за кожей лица</b><br>Инструкция по применению   | 61-71 |
| <b>(PL) Urządzenie do masażu<br/>ultradźwiękami</b><br>Instrukcja obsługi | 72-82 |

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## Dear customer

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of beauty, heat, gentle therapy, blood pressure/diagnosis, weight, massage and air.

Please read these instructions for use carefully and observe the information they contain.

With kind regards,  
Your Beurer team

|                                      |    |
|--------------------------------------|----|
| 1. Intended use.....                 | 12 |
| 2. Safety notes .....                | 13 |
| 3. Getting to know the device.....   | 14 |
| 4. Device description .....          | 15 |
| 5. Initial use .....                 | 16 |
| 6. Operation.....                    | 16 |
| 7. Cleaning and maintenance.....     | 19 |
| 8. What if there are problems? ..... | 20 |
| 9. Technical data .....              | 20 |
| 10. Disposal.....                    | 21 |

### Warning!

Non-observance can lead to permanent injury.

### Caution!

Non-observance can lead to minor injuries or damage.

### Important!

Non-observance can lead to material or environmental damage.

Cross-reference to another section.

## Included in delivery

- Ultrasound device with cover
- Charger
- Power supply unit with mains cable
- Cream (15ml) (You can reorder the moisturising cream specially designed for use with the FC80 / „Pureo Sonic“ from your local retailer using order number 162.944. To do so, please note the service address list supplied.)

## 1. Intended use

The ultrasound device is intended for the cosmetic treatment of facial skin. Please observe the safety notes in these instructions for use.

Use the device only according to the methods specified.

The device is intended for private use.

## 2. Safety notes

- Please read the instructions for use carefully! Non-observance of the following information may result in personal injury or material damage.
- Store these instructions for use in a safe place and make them accessible to all users.
- Make sure you include these instructions for use when handing over the device to third parties.

### Electric shock

#### Warning!

- As with every electrical device, use this ultrasound device carefully and cautiously in order to avoid hazards due to electric shock.
- Use device only with the power supply unit provided and at the mains voltage specified on the type plate.
- Only plug the power supply unit into the socket when the device is switched off.
- Never use the device if it or its accessories are visibly damaged.
- Do not use the device if it has sustained a fall or if water has penetrated the device interior.
- Do not use the device during a storm.
- In the event of defects or operational faults, switch off the device immediately. Pull the plug out of the mains socket.
- Do not pull on the mains cable or on the device to disconnect the plug from the socket.
- Never hold or carry the device by the mains cable.
- Do not clamp or pinch the mains cable.
- Keep the mains cable away from hot surfaces.
- Ensure that the mains cable and plug do not come into contact with water, steam or other liquids.
- Only use the device in dry internal rooms (e.g. never in the bathroom, sauna, etc.).
- Never reach for a device that has been submerged in water. Unplug the power supply unit from the socket immediately.

- Do not subject the device to any impacts or drop the device.

### Repairs

#### Warning!

- Do not open the device.
- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact customer service or an authorised retailer.
- If the power supply unit or mains cable of the device is damaged, it must be replaced by the manufacturer or an authorised customer service point.

### Risk of fire

#### Warning!

There may be a risk of fire if the device is used improperly or these instructions for use are ignored!

- Never use the device underneath a cover, such as a blanket or pillow.
- Never use the device near highly flammable substances.

### Operation

#### Warning!

- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.
- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or are instructed by him in how to use the device.

- Children should be supervised around the device to ensure they do not play with it.
- The device must be switched off at the On/Off switch after every use and before cleaning. Unplug the power supply unit from the mains socket.
- Ensure that the mains cable does not pose a trip hazard.
- Do not treat the following areas: the eye and its immediate vicinity, and the ear.
- Use the ultrasound device only on healthy and intact skin.
- Do not use the ultrasound device if you have electrical or metal implants.

### **Caution!**

- Do not place any objects on the device.
- Protect the device against high temperatures.


## **3. Getting to know the device**

With the Beurer FCE80 you can perform a cosmetic facial treatment yourself at home. What was previously only possible in the cosmetic studio, you can now do at home on a regular basis using the FCE80.

The ultrasound head generates around 5 million gentle oscillations per second. These oscillations penetrate into your skin, thus stimulating the various layers of skin. The movement reaches right down to the deep skin layers, makes these

more permeable and transports any care products used and moisture down to these deep areas. This stimulates the metabolism and aids detoxification. The ultrasound treatment has an effect upon the lymphatic system. Your cells regenerate more quickly. Taken together, these factors can improve the overall appearance of your skin, leaving you looking younger and healthier.

## 4. Device description

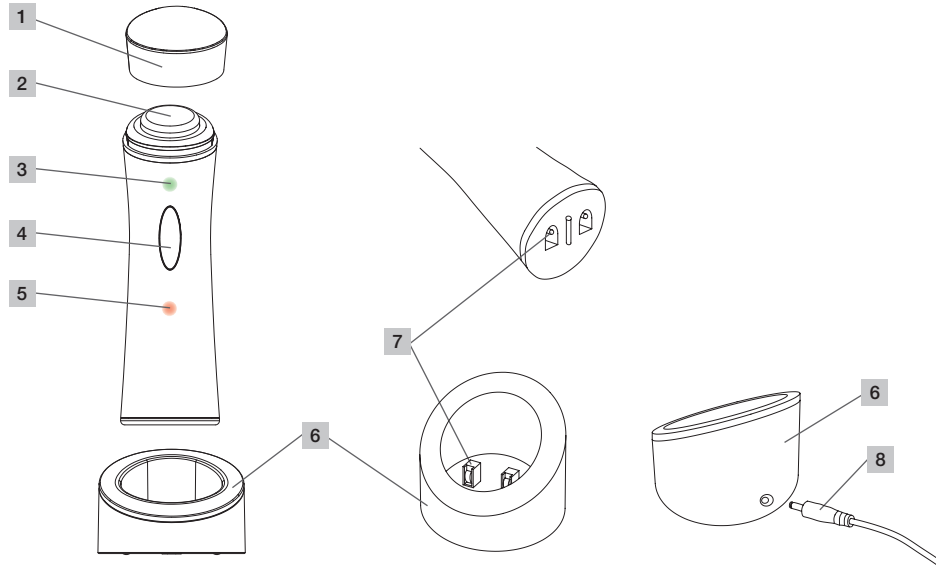1 Transparent cover

2 Ultrasound head

3 Green LED,  
lights up when the device is in use

4 On/Off switch

5 Red LED,  
lights up while the battery is charging

6 Charger

7 Battery contacts

8 Charger plug

## 5. Initial use

### Unpacking the device

---

#### **Warning!**

Keep packaging material away from children!  
There is a choking hazard.

---

- > Open the cardboard packaging
- > Remove the ultrasound device with transparent cover, the charger and the cable with the power supply unit
- > Check all parts for damage

### Installation location

---

#### **Warning!**

- Before using the device for the first time, read the safety notes in the Safety chapter.
    - ▷ Safety pages 13 to 14
  - Ensure that the mains cable does not pose a trip hazard.
- 

- > Insert the mains cable plug into the charger
- > Place the charger on a smooth, level surface
- > Plug the power supply unit into a mains socket

## 6. Operation

---

#### **Warning!**

- Before using the device for the first time, read the safety notes in the Safety chapter.
    - ▷ Safety pages 13 to 14
- 

#### **Caution!**

- If you experience irritations of any kind, stop treatment immediately.
  - When using the ultrasound device for the first time, or if you have sensitive skin, reduce the treatment time.
- 

### Charging the battery

*Place the device in the charger*

*The red LED lights up; this means that the charging process has begun and the battery is being charged.*

The entire charging process takes around 3 hours for the first charge or whenever the device has not been used for a relatively long time. When the device has been in regular use, charging takes around 2 hours.

*The red LED goes out as soon as the charging process has finished.*

---

#### **Important!**

To ensure the long life of your battery, the device should always be fully discharged before being fully recharged (takes around 2 hours).

---

### Checking that the ultrasound is working

Check that the ultrasound head is working.

- > Carefully let a few drops of water drip onto the ultrasound head
- > Switch on the device using the On/Off switch  
*If the drops of water vibrate, the device is working.*

#### **Important!**

- Clean your face thoroughly before treatment.
  - Use only water-based gels and skincare products. Oily skincare products do not support the function of the ultrasound device.
  - Check whether your skin tolerates the skincare product/contact gel you are using.
-

## Skin compatibility test

Simply test the compatibility of the skincare product so that you can be sure you will have no skin problems.

- > Clean a patch of skin in the crook of your arm thoroughly with soap, apply a little of the skincare product to this cleaned patch of skin and allow the product to soak in and dry.

- If the skin shows no reactions within the next 48 hours, you can use the tested product without hesitation.
- If skin reactions appear before the end of the 48-hour period, you should not use the tested product for facial treatment.

## Risks and side effects

---

### **Warning!**

- The device should not be used by people with heart conditions, active vascular diseases, pacemakers, thrombosis, liver or kidney diseases, serious lipid metabolic diseases or by people taking blood thinning medications.
  - Use is not recommended for people who have received medical transplants, individuals who are pregnant or breastfeeding, or people with metal prostheses or terminal diseases.
- 

---

### **Important!**

- The ultrasound treatment stimulates circulation and flushes out waste products. This cleansing effect can, in some cases, result in the skin feeling less firm. This is particularly applicable for skin types with thin connective tissue. If this occurs, treatments should be performed at intervals of at least one week.
- 
- The use of the care concentrate may trigger an allergic reaction.
- 

In case of uncertainty, always consult your doctor before use.

## Recommended treatment period and frequency

We recommend facial treatment a maximum of 3 times per week for 10 minutes per treatment.

## Beginning the treatment

- > Remove the transparent cover from the ultrasound head
- > Distribute some skincare product over your face
- > Press the On/Off switch to switch the device on  
*The green LED lights up*

## Moving the ultrasound head


The movements of the ultrasound head over the face should correspond approximately to the natural lymphatic flow. Follow one of the recommended sequences for facial treatment.

- Put the ultrasound head on your skin. Press the device down gently while you let it glide over your face.
- Guide the device slowly from the centre of your face outwards and from the bottom of your face to the top.


- Lift the device as soon as you arrive at the end of a movement and begin again at the starting point.
- Always complete a cycle.
- Do not move the device back and forth.
- Move the device constantly during treatment; do not hold it in one position for long periods.

## Recommended sequence for full facial treatment


- Duration of each movement approx. 5 seconds.
- Repeat each movement eight times.
- Complete a cycle, then perform the entire cycle for the other half of the face.

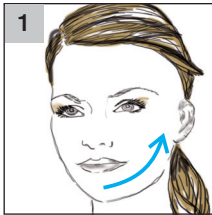
- 1 > Move from your chin to below your ear

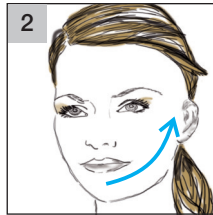
- 2 > Move from your chin to in front of your ear

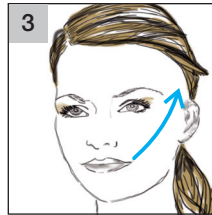
- 3 > Move from the corner of your mouth to above your ear

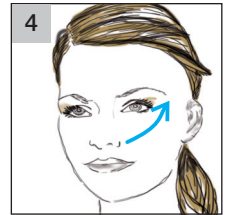
- 4 > Move from below your nose to above your ear

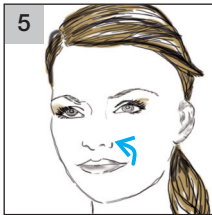- 5 > Move from in front of the corner of your mouth to your nostril

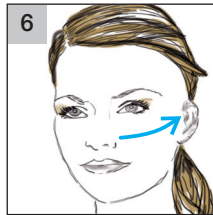- 6 > Move from your nostril to your cheekbone

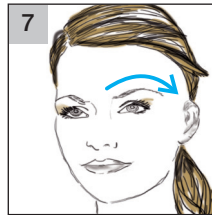- 7 > Move from the centre of your forehead to the side in an arc

## 7. Cleaning and maintenance

---

### **Warning!**

- Observe the safety notes in the Safety chapter.
  - ▷ Safety pages 13 to 14
- The device must be switched off at the On/Off switch after every use and before cleaning.
- Unplug the power supply unit from the mains socket.
- Never clean the charger with the power supply unit plugged in.

- 
- > Switch the device off using the On/Off switch
  - > Unplug the power supply unit from the mains socket
  - > Wipe the ultrasound head with a tissue  
*Skincare product residues should be thoroughly removed.*

- > If necessary, clean with a slightly damp cloth

---

### **Important!**

- Clean the device using only the methods specified.
- Under no circumstances may liquid enter the charger.
- Do not use any solvent-based cleaning products.
- Only use a mild flushing agent or vinegar/vinegar-based cleaner to perform cleaning or disinfection work.

---

### **After cleaning**

- > Plug in the power supply unit
- > Charge the device
- ▷ Initial use ▷ Charging the battery

## 8. What if there are problems?

| Problem  | Cause   | Solution   |
|--|---|--|
| The device does not respond after switching on.  | Is the battery flat?  | <ul style="list-style-type: none"> <li>&gt; Charge the battery</li> <li>▷ Initial use</li> <li>▷ Charging the battery</li> </ul> |
| The green LED is flashing.   | Is the battery almost flat?   | <ul style="list-style-type: none"> <li>&gt; Charge the battery</li> <li>▷ Initial use</li> <li>▷ Charging the battery</li> </ul> |
| Fault during battery charging. After putting the device on the charger, the red LED does not light up. | Is the power supply unit not plugged in correctly?  | > Plug the power supply unit into the mains socket   |
|  | Is the plug on the mains cable not correctly inserted into the charger?                                     | > Insert the plug into the back of the charger   |
|  | Is the device not correctly positioned on the charger?  | > Remove the device and replace it on the charger. The charger must be standing up straight and the device may not tip over.     |
|  | Are the battery contacts of the device and charger dirty?   | > Using a lint-free cloth, clean the contacts of the charger and device  |
| The device is not fully charged after 2 hours. The red LED is lit.                                     | Is the device being used for the first time? Has it been six months or more since the device was last used? | > If the battery has been particularly heavily discharged, charging can take longer than 2 hours                                 |

## 9. Technical data

|                           |  |
|---------------------------|--|
| Dimensions (diameter x h) | Approx. 7 cm x 17.5 cm with charger  |
| Mains voltage/frequency   | 100-240 V / 50/60 Hz   |
| Battery                   | NiMH   |
| Battery life              | Approx. 20 minutes per charge, life cycle approx. 500 charges  |
| Signs and symbols         |  The unit is double-protected and thus corresponds to protection class 2. |

## 10. Disposal

The rechargeable NiMH battery is a valuable resource and must be disposed of separately.

### Removing the battery before disposal

- > Switch the device off using the On/Off switch  
*The green LED must not light up*
- > Remove the screw on the base of the device
- > Remove the base
- > Withdraw the ultrasound head with battery holder from the housing

- > Withdraw the battery plug from the circuit board
- > Remove the NiMH battery and dispose of it separately

Observe the local regulations for material disposal. Dispose of the unit in accordance with EC Directive 2002/96/EC – WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal.

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