(D) Warmwachs-Haarentferner Gebrauchsanleitung	(RUS) Прибор для удаления волос воском Инструкция по применению
(GB) Warm wax hair remover Instruction for Use	(PL) Depilator woskowy Instrukcja obsługi
(F) Épilateur à la cire chaude Mode d'emploi.....	(NL) Harsapparaat Gebruikshandleiding.....
(E) Aplicador de cera depilatoria tibia Instrucciones para el uso.....	(P) Depilador a cera quente Instruções de utilização
(I) Depilatore a cera calda Istruzioni per l'uso	(GR) Αποτριχωτική μηχανή με ζεστό κερί Οδηγίες χρήσης.....
(TR) Sıcak Ağdalı Epilasyon Cihazı Kullanma Talimatı	

Inhalt

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ENGLISH

Contents

1. Proper use.....	8	5. Care and storage.....	11
2. Safety information	8	6. Replacement parts and wearing parts....	11
3. Appliance description	10		
4. Operation.....	10		

Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications related to heating, weight, blood pressure, body temperature, pulse, gentle therapy, massage, air and beauty.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,
Your Beurer Team

WARNING

- The unit is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risks.

Items included in the package:

- Base station
- Base device/applicator
- 2 wax cartridges
- 50 fabric strips
- These instructions for use

Explanation of symbols

The following symbols are used in the operating instructions:

WARNING

Warning of risks of injury or health hazards

CAUTION

Safety information about possible damage to appliance/accessories

NOTE

Important information

CAUTION

This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – danger of electric shock!

The device has double protective insulation and therefore complies with protection class 2

Only use indoors.

1. Proper use

The warm wax hair remover enables you to carry out thorough and quick hair removal in your own home, while simultaneously adding moisture to your skin and making it softer.

The applicator can be used to apply the wax to the skin quickly and easily, in thin layers, without burning. The base station enables two wax cartridges to be warmed simultaneously to ensure continuous hair removal without waiting periods.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

This device is intended for household use only, not for commercial applications.

2. Safety information

Please read these operating instructions carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the appliance if you sell it or pass it on.

WARNING

- Keep children away from the packaging materials. Risk of suffocation.
- Children must be supervised when using this device.

CAUTION

- Before use, make sure that the device and accessories do not show any visible signs of damage. If you are in any doubt, do not use it and contact your dealer or a specified customer service address.
- Do not stretch, twist or kink the power cord.
- Do not position the power cord on or pull it over sharp or pointed objects.
- Only insert the plug into a socket that has been installed in compliance with applicable regulations.
- Do not open the case under any circumstances.
- Keep the device away from pointed or sharp objects.

- If the device has been dropped or otherwise damaged, it must not be used again.
- Switch the device off and unplug it from the power supply after every use and before every cleaning.

WARNING: Electric shock

As with every electrical device, use this device carefully and cautiously, in order to avoid hazards due to electric shock.

For this reason, the device must only be operated as follows:

- only use with the mains voltage noted on the base station.
- never use the device if any damage to the appliance itself or to any of its accessories is evident
- never during a thunderstorm.

In the case of defects or operational faults, switch the device off immediately and disconnect the device from the mains socket. Do not pull on the mains cable or on the device to disconnect the plug from the socket. Never hold or carry the device by the mains cable. Keep the cables away from hot surfaces.

Keep the device, the plug and cables away from water or other liquids.

For this reason, the device must be operated as follows:

- only in dry inside rooms
- only with dry hands
- never in the bath, in the shower or near a filled basin
- never in swimming pools, whirlpools or in the sauna
- never out in the open

Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the mains outlet.

WARNING : Repairs

- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorised dealer.
- In the event that the mains cable is damaged, this must be replaced by an authorised retailer.
- To order replacement parts (wax cartridges and fabric strips), please refer to the enclosed service address list.

WARNING: Fire hazard

There is a risk of fire if the device is used incorrectly or if it is not used according to these operating instructions.

For this reason, the massage device must only be operated as follows:

- Never leave the device unsupervised, particularly if children are around.
- Never use the device under covers, such as blankets or pillows.
- Never use the device in the vicinity of petrol or other easily flammable materials.


CAUTION: Disposal

Follow local regulations regarding the disposal of materials. Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please contact the appropriate local authorities.

3. Appliance description

1. Base station
2. Charger contacts
3. Base station protective cap
4. Base device/applicator
5. Indicator lamp
6. Fill level display
7. Protective cap for applicator
8. Warm wax cartridges (2)
9. On/Off switch

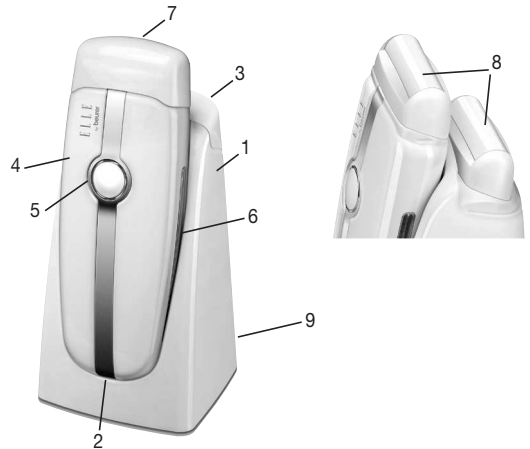4. Operation

- Remove the packaging.
- Check the device, power supply and cables for damage.
- Connect the device to the mains power. The device must be switched off at this time.
- Route the cable safely so that there is no risk of anyone tripping over it.
- Switch on the device using the On/Off switch.
- Place the applicator on the base station and insert the second wax cartridge into the separate heating chamber as shown in the image.
- The indicator lamp flashes while heating is taking place. Once the wax has heated fully, the indicator lamp then remains continuously lit. This process lasts approx. 25–30 minutes.
- As soon as the wax is hot, roll the wax roller until it rotates without any problems.
- If you are using a warm wax hair remover for the first time, please refer to the instructions to ensure you are using it correctly.
- In order to become accustomed to the hair removal system and to ensure that you are not allergic to it, it is preferable to start with a small area on the lower part of the leg. Then wait 24 hours before carrying out the complete hair removal process.
- The skin must be clean and dry prior to carrying out hair removal.

NOTE

- The wax cartridge contains 50 ml wax. This quantity is sufficient for two half legs and a small number of follow-up treatments. The more practice you get the less wax you will need to use.
 - In order to ensure optimum hygiene and cleanliness, the wax roller and the warm wax cartridge have been designed for single use and must be disposed of following use. Doing so prevents infections and inflammations.
 - Never remove hair straight after a bath or shower; wait approx. two hours to ensure the skin is not irritated.
 - If you want to use a sunbed within 24 hours of carrying out hair removal, protect your skin by using sunscreen.
-
- Apply the wax in the direction of hair growth, ensuring that the roller rotates and that the device is held at a 45° angle to prevent wax leaking out.
 - Insert the applicator onto the base station between uses to ensure that the wax retains its temperature and does not leak out.
 - The fabric strip must be laid on the skin immediately after the wax is applied. Lay the strip in the direction of hair growth, covering the entire area, to ensure that the strip sticks to the wax effectively.
 - Keep the skin taut using your free hand and pull away the strip using a firm, backward movement, against the direction of hair growth and parallel to the surface of the skin. Do not pull diagonally or upwards as this may cause the hair to break and not be pulled out.

- Each paper strip can be used several times. The more often a strip has been used, the more firmly it sticks.
- If the first cartridge has been used up, continue the hair removal process using the replacement cartridge. To do so, remove the cartridge from the chamber with caution, in order to avoid getting burnt, and insert this into the applicator.

CAUTION

- Only use the warm wax hair remover following the method described in this manual. Only use beurer wax cartridges.
- Do not treat any sensitive areas such as eyebrows, birthmarks, warts, open wounds, cuts, burns etc.
- This device is not suitable for persons who suffer from diabetes, varicose veins or circulation problems, or persons who are prone to infections.
- Do not use the wax on irritated skin. Do not apply the wax to areas to which other substances have already been applied.
- Please do not use any wax on the face area, ear lobes, the inside of the ears, on nipples, or on irritated, cracked or sunburnt skin.
- If any wax remains on the skin following hair removal, the strip must be reapplied, rubbed down firmly and pulled off. The remaining wax can also be carefully removed using oil wipes.
- If you want to use the warm wax hair remover again, heat the wax until it is liquid. Hold the applicator downwards until the remaining wax is hot on the roller. If the wax does not completely liquefy, roll the wax roller along a fabric strip until the wax roller rotates completely.
- Once you have completed the treatment, place the base device on the base station and leave the wax to cool. Always position the device upright, to ensure that the wax does not leak out.

5. Care and storage

WARNING

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Do not use the device again until it is completely dry.
- Do not clean the device in the dishwasher.
- Do not use aggressive cleaning agents or hard brushes.

First unplug the device and allow it to cool off. Then you can clean it with a soft, slightly moistened cloth. Then wipe it dry with a soft, lint-free cloth.

6. Replacement parts and wearing parts

	Item number and/or order number
2x wax cartridges + 50 fabric strips	576.41*

*Commercially available or from the corresponding service address (according to the service address list)

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